
CORIO CD-201F Refrigerated / heating circulator


Refrigerated / heating circulators of the new CORIO series distinguish themselves with a great price-to-performance ratio. They are ideal for all standard tasks and routine work in laboratories and industry.

Product features

- Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- Easy pump change-over between internal and external circulation
- External pump connections
- USB connection
- Space-saving cooling coil design provides more usable space in the bath tank
- Bath lid and drain tap included
- Removable ventilation grid
- Refrigeration unit without side vents
- Class III (FL) according to DIN 12876-1
- Very quiet
- High-quality bath tanks made of stainless steel with bath lid and drain tap

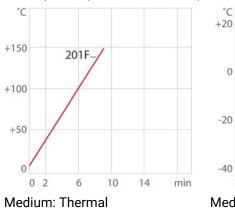

Heat-up time (230V, 2kW heater) Cool-down time

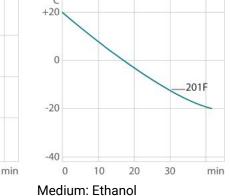

Performance values

115V/60Hz (Nema N5-15 Plug)	
Heating capacity kW	1
Viscosity max. cST	50
Pump capacity flow rate l/min	16
Pump capacity flow pressure psi	4.8
Power A	12



Julabo



Order No.			9012702.02	
Cooling capacity (Ethan	iol)			
°C	20		0	-20
kW	0.22		0.16	0.06
				with ethanol; over 20 °C with thermal oil e values may differ with other bath fluids.
Refrigerant stage 1				
Refrigerant	R134a			
Filling volume g	70			
Global Warming Potenti for R134a	ial 1430			
Carbon dioxide equivale t	ent 0.1			
230V/60Hz (Schuko Plug - CE	EE 7/4 Plug Type F)			
leating capacity kW				2
/iscosity max. cST				50
Pump capacity flow rate I/min				17
Pump capacity flow pressure ps	i			6.2
Power A				12
Order No.			9012702.13	
Cooling capacity (Ethan	iol)			
°C	20	0	-10	-20
kW	0.22	0.16	0.12	0.06
				with ethanol; over 20 °C with thermal oil a values may differ with other bath fluids.
Refrigerant stage 1				
Refrigerant	R134a			
Filling volume g	70			
Global Warming Potenti for R134a	ial 1430			
Carbon dioxide equivale	ont 01			

Carbon dioxide equivalent 0.1 t

Technical data

Available voltage ve	ersions	Bath	
Order No.	9 012 702	Bath tank	Stainless steel
Available voltage version	ns:	Bath cover	integrated
9012702.01	100V/50-60Hz (Nema N5-15 Plug) (R134a)	Usable bath opening in. (W x L / D) $$	5.1 x 5.9 / 5.9
9012702.02	115V/60Hz (Nema N5-15 Plug) (R134a)		
9012702.03	230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F) (R134a)		



9012702.03.chn	230V/50Hz (CN Plug) (R134a)
9012702.04	230V/50Hz (UK Plug Type BS1363A) (R134a)
9012702.05	230V/50Hz (CH Plug Type SEV 1011) (R134a)
9012702.13	230V/60Hz (Schuko Plug - CEE 7/4 Plug Type F) (R134a)
	(R134a) 230V/60Hz (Schuko Plug - CEE 7/4 Plug

Cooling		Other	
Cooling of compressor	1-stage Air	Classification	Classification III (FL)
	IP Code	IP 21	
		Pump function	Pressure Pump
	Pump type	Immersion Pump	

Electronics	
Digital interface	USB
Temperature control	PID1
Absolute temperature calibration	1 Point Calibration
Temperature displayTemperature display	LED
Temperature settingTemperature setting	Keypad

	•
Dimensions and volumes	
Weight lbs	55.1
Barbed fittings inner diameter	8/12 mm
Dimensions in. (W \times L \times H)	17.3 x 16.1 x 17.3
Filling volume I	3 4
Pump connections	M16x1 male

Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male) $\,$

All Benefits

display °C

Temperature values

Temperature stability °C

Ambient temperature °C

Working temperature range °C

Temperature display resolution °C

Setting the resolution of the temperature



100% Checked.

100% testing. 100% quality. Each JULABO Circulator undergoes thorough quality testing before leaving the factory.

0.1

-20 ... +150 ±0.03

+5 ... +40

0.01 ... 0.1



Handle with ease.

Makes day-to-day work easy. Comfortably move your CORIO around by using the ergonomic handles (front and rear).



Space saving. Free up space.

Place your JULABO Circulator right next to an application, another unit, or wall. That saves space. This is made possible by eliminating vents and connections on the sides.



ATC. Absolute Temperature Calibration, 1-point calibration (CD).

More bath.

Designed for more comfort. Thanks to the recessed cooling coil, the internal bath provides more space.



Condensation protection. Superb design solution. Integrated ventilation directs air over the bath lid and minimizes

directs air over the bath lid and minimizes condensation.





Green technology.

Development consistently applied environmentally friendly materials and technologies.



Modern. Reliable. High-grade components in every JULABO

Circulator - platinum sensors, proven motor technology, CAN-Bus communication and much more.



Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



Solid. Minimized energy loss through high-quality insulation.



Tidy. The special drain tap for easy draining of bath fluids without tools.



JULABO. Quality. Highest standards of quality for a long product life.



Quick start. Individual JULABO consultation and comprehensive manuals at your disposal.



Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes.



Connection. Easy.

Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.

Exact.

You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO precise and perfect.



Switch on. And off you go. Intelligent operating concept. Ready for operation

with just a few quick and easy steps.

Everything at the front. All operating controls and safety functions are accessed easily and comfortably from the front.



Locked in. The lockable power plug guarantees safe connection. More process safety.



Safety. CORIO CD and CP comply with Class III (FL) according to DIN 12876-1 and switches off automatically in case of high temperature or low liquid level alarm.

The pump is controlled via a lever located directly

below the display. Easily change between internal

Services 24/7.

Internal. External.

and external circulation.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.us.


Stable.

Rubber feet allow for a secured footing of your CORIO to prevent damage to your laboratory equipment.




Touching permitted. Optimum safety. The ergonomic plastic handle

protects your fingers from hot surfaces.

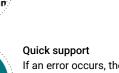

Brilliant.

a distance.

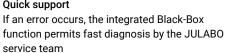

Very bright display makes it easy to read even from







Connectivity.

Remote control made easy. CORIO CD circulators feature a USB connection.



Early warning system for low liquid level Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down