


CORIO CD-B19 Heating circulator with stainless steel open bath

The open heating bath circulators of the new CORIO series feature high-quality bath tanks made of stainless steel.

Product features

- · Models for internal applications
- Bright, white, easy to read display
- Very quiet
- Easy switching between internal and external circulation
- USB connection
- High-quality bath tanks made of stainless steel
- Integrated drain screw (except B5)
- Integrated drain screw

Heat-up time (230V, 2kW heater)

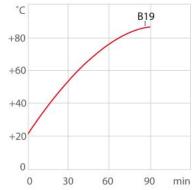

Medium: Water

Performance values

| 115V/60Hz (Nema N5-15 Plug) | |
|---------------------------------|-----|
| Heating capacity kW | 1 |
| Viscosity max. cST | 50 |
| Pump capacity flow rate I/min | 15 |
| Pump capacity flow pressure psi | 4.8 |
| Power A | 10 |

| 230V/60Hz (Schuko Plug - CEE 7/4 Plug Type F) | |
|---|-----|
| Heating capacity kW | 2 |
| Viscosity max. cST | 50 |
| Pump capacity flow rate I/min | 17 |
| Pump capacity flow pressure psi | 6.2 |
| Power A | 10 |





Technical data

| Available voltage ver | rsions | Bath | |
|----------------------------|--|-------------------------------------|-------------------|
| Order No. | 9 012 419 | Bath tank | Stainless steel |
| Available voltage versions | s: | Bath cover | accessory |
| 9012419.01 | 100V/50-60Hz (Nema N5-15 Plug) | Usable bath opening in. (W x L / D) | 11.8 x 13.8 / 5.9 |
| 9012419.02 | 115V/60Hz (Nema N5-15 Plug) | | |
| 9012419.03 | 230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F) | | |
| 9012419.03.chn | 230V/50Hz (CN Plug) | | |
| 9012419.04 | 230V/50Hz (UK Plug Type BS1363A) | | |
| 9012419.05 | 230V/50Hz (CH Plug Type SEV 1011) | | |
| 9012419.13 | 230V/60Hz (Schuko Plug - CEE 7/4 Plug Type F) | | |

| Cooling | | Other | |
|------------------------|----------------|-------------------------|--|
| Cooling coil accessory | Classification | Classification III (FL) | |
| | IP Code | IP 21 | |
| | Pump function | Pressure Pump | |
| | Pump type | Immersion Pump | |

| Electronics | | Dimensions and volumes | |
|--|---------------------|--|------------------|
| Digital interface | USB | Weight lbs | 25.4 |
| Temperature control | PID1 | Barbed fittings inner diameter | 8/12 mm |
| Absolute temperature calibration | 1 Point Calibration | Dimensions in. $(W \times L \times H)$ | 15 x 22.8 x 16.5 |
| Temperature displayTemperature display | LED | Filling volume I | 14 19 |
| Temperature settingTemperature setting | Keypad | Pump connections | M16x1 male |

| Temperature values | | Included in delivery |
|--|----------|--|
| Setting the resolution of the temperature display °C | 0.1 | 2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male) |
| Working temperature range °C | +20 +150 | |
| Temperature stability °C | ±0.03 | |
| Ambient temperature °C | +5 +40 | |
| Temperature display resolution °C | 0.01 0.1 | |

 $For applications \ near \ or \ below \ ambient \ temperature: use \ a \ cooling \ coil \ or \ JULABO \ immersion \ cooler.$

All Benefits



100% Checked.

100% testing. 100% quality. Each JULABO Circulator undergoes thorough quality testing before leaving the factory.



ATC.

Absolute Temperature Calibration, 1-point calibration (CD).



Green technology.

Development consistently applied environmentally friendly materials and technologies.



Internal. External.

The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.





JULABO. Quality.

Highest standards of quality for a long product life.



Modern. Reliable.

High-grade components in every JULABO Circulator – platinum sensors, proven motor technology, CAN-Bus communication and much more.



Satisfied customers

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.us.



Stable.

Rubber feet allow for a secured footing of your CORIO to prevent damage to your laboratory equipment.



Quick start.

Individual JULABO consultation and comprehensive manuals at your disposal.



Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes.



Brilliant.

Very bright display makes it easy to read even from a distance.



Connection. Easy.

Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



Everything at the front.

All operating controls and safety functions are accessed easily and comfortably from the front.



Exact

You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO precise and perfect.



Locked in

The lockable power plug guarantees safe connection. More process safety.



Switch on. And off you go.

Intelligent operating concept. Ready for operation with just a few quick and easy steps.



Early warning system for low liquid level

Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



Connectivity.

Remote control made easy. CORIO CD circulators feature a USB connection.