




Product data sheet

VALEGRO 1201 Recirculating cooler

Modern, user-friendly chillers for a wide range of cooling tasks

The air-cooled VALEGRO recirculating coolers are environmentally friendly all-rounders for a wide range of cooling and heating tasks from -20 to +40 °C. The pump can be user-friendly adjusted in different pump stages on a percentage basis.

Whether in the laboratory, research, or production, the VALEGRO 1201 recirculating coolers offer a constant cooling capacity of 1.2 kW at 20 °C, supporting stable processes around the clock. These modern devices offer precise temperature control for a wide range of applications and, thanks to their compact design, fit seamlessly into any working environment. The refrigeration circuit uses an eco-friendly, natural refrigerant. When designing the VALEGRO recirculating coolers, the focus was on maximum user-friendliness. Many helpful features such as a timer function, a removable venting grid that can be removed without tools, and an illuminated and simple to clean level indicator make daily work easier. The inclined OLED touch display is easy to read and ergonomic to use. The intuitive menu navigation allows quick access to all functions. Warning notices and system messages are displayed in clear, understandable language, ensuring maximum safety when using the device. In addition to RS232 and USB-C interfaces, VALEGRO recirculating coolers are also available with Ethernet (optional). All interfaces are located directly below the display and are easy to use from the front.

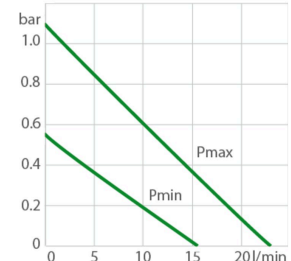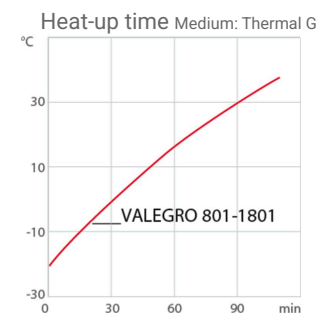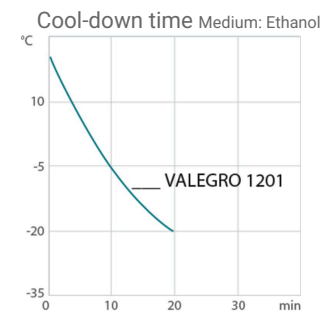


Product Features

- ✓ User-friendly operation thanks to comprehensive help texts
- ✓ Intuitive menu navigation in full text German/English
- ✓ Data recording
- ✓ Autostart function
- ✓ Simple timer function
- ✓ Removable venting grid without tools
- ✓ Precise temperature stability of ± 0.3 °C
- ✓ Casters for easy handling

Available options

Interfaces

Ethernet
VALEGRO Ethernet activation



Performance values

| | |
|--------------------------------------|------------|
| 115V/60Hz (Multiple Plugs (EU / US)) | |
| Viscosity max. cSt | 50 |
| Pump capacity flow pressure l/min | 15 ... 23 |
| Pump capacity pressure psi | 7.3 ... 16 |
| Power consumption A | 7 |
| Voltage Tolerance % | 115V ±10 |

Order number including voltage version (115V/60Hz) 9610121.30005627

Cooling capacity 1 (Ethanol)

| | | | | |
|-----------------|-----|-----|-----|------|
| °C | 20 | 0 | -10 | -20 |
| kW ¹ | 1.2 | 0.9 | 0.6 | 0.35 |

Note about natural refrigerants:
Temperature control units using natural refrigerants are often subject to regulatory requirements regarding the installation site, operation, transport or disposal of the units. If you have any questions, we will be happy to advise you.

Refrigerant stage 1

| | |
|-----------------------------------|---------|
| Refrigerant | R290 |
| Filling weight g | 80 |
| Global Warming Potential for R290 | 3 |
| Carbon dioxide equivalent t | 0.00024 |

¹ Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

Technical data

Available voltage versions

Order No. 9610121

Available voltage versions:

9610121.30005627 100-230V 50/60Hz (Multiple Plugs (EU / US)) (R290)

9610121.30005628 100-230V 50/60Hz (UK Plug Type BS1363A) (R290)

9610121.30005629 100-230V 50/60Hz (Multiple Plugs (EU / CH)) (R290)

9610121.30005630 100-230V 50/60Hz (CN Plug) (R290)

Cooling

Cooling of compressor 1-stage Air

Bath

Bath tank Stainless steel

Electronics

Interfaces Ethernet optional, RS232, USB

Temperature control PID1

Temperature display OLED

Temperature setting Touchscreen

Electronic Timer h:min 99 ... 59

Dimensions and volumes

Weight lbs 94.8

Barbed fittings inner diameter mm 8/12 mm

Total dimensions in. (W × L × H) 14.4 x 20.7 x 24

Body dimensions in. (W × L × H) 14.4 x 18.7 x 23.6

Filling volume l 4 ... 7

Pump connections M16x1 male

Temperature values

Return flow temperature max. °C +85

Working temperature range °C -20 ... +40

Temperature stability °C ±0.3

Ambient temperature °C +5 ... +40

Temperature display resolution °C 0.1

Other

Classification Classification I (NFL)

Standards CE, RoHs

Alarm notification Acoustic, Optical

IP Code IP 21

Pump function Pressure Pump

Pump type Centrifugal Pump

User Interface Language English, German

Sound pressure level dbA² 42 ... 51

² The sound pressure level may increase as the ambient temperature rises.

Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections
M16x1 male)

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio. Technical changes without prior notification. Images may deviate from the original. | Datasheet No.EN-US9610121/251125