

Product data sheet

VALEGRO 1803 recirculating cooler

with cooling capacities of 1.2 - 2.5 kW

The more powerful VALEGRO recirculating coolers are designed for larger and more complex applications. They feature a precisely adjustable, increased pump pressure up to 3.5 bar.


VALEGRO 1803 Whether in laboratories, research, or production, VALEGRO recirculating coolers offer a constant cooling capacity of 1.8 kW at 20 °C, supporting stable processes around the clock. These modern units offer precise temperature control for a wide range of applications and, thanks to their compact design, fit seamlessly into any working environment. The refrigeration circuit uses an eco-friendly, natural refrigerant. When designing the VALEGRO recirculating cooler, the focus was on maximum user-friendliness. Many helpful features such as a timer function, a removable venting grid without tools, and an illuminated and simple to clean level indicator make daily work easier. The inclined OLED touch display is easy to read and ergonomic to use. The intuitive menu navigation allows quick access to all functions. Warning notices and system messages are displayed in clear, understandable language, ensuring maximum safety when using the device. In addition to RS232 and USB-C interfaces, VALEGRO recirculating coolers are also available with Ethernet (optional). All interfaces are located directly below the display and are easy to use from the front.



Cool-down time Medium: Ethanol



Heat-up time Medium: Thermal G

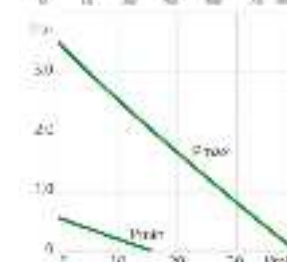

Product Features

- ✓ Adjustable pump pressure
- ✓ User-friendly operation thanks to comprehensive help texts
- ✓ Intuitive menu navigation in full text German/English
- ✓ Data recording
- ✓ Autostart function
- ✓ Simple timer function
- ✓ Removable venting grid without tools
- ✓ Precise temperature stability of $\pm 0.3\text{ }^{\circ}\text{C}$

Available options

Interfaces

Ethernet
VALEGRO Ethernet activation



Performance values

230V/60Hz (Multiple Plugs (EU / US))

Viscosity max. cSt	50
Pump capacity flow pressure l/min	15 ... 40
Pump capacity pressure psi	7.3 ... 50.8
Power consumption A	7
Voltage Tolerance %	230V ±10

Order number including voltage version (230V/60Hz) 9610183.30006991

Cooling capacity 1 (Ethanol)

°C	20	0	-10	-20
kW ¹	1.8	1.3	0.8	0.4

Note about natural refrigerants:

Temperature control units using natural refrigerants are often subject to regulatory requirements regarding the installation site, operation, transport or disposal of the units. If you have any questions, we will be happy to advise you.

Refrigerant stage 1

Refrigerant R290

Filling weight g 99

Global Warming Potential for R290 3

Carbon dioxide equivalent t 0.000297

¹ Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

Technical data

Available voltage versions

Order No. 9610183

Available voltage versions:

9610183.30006991 200-230V 50/60Hz (Multiple Plugs (EU / US)) (R290)

9610183.30006992 200-230V 50/60Hz (UK Plug Type BS1363A) (R290)

9610183.30006993 200-230V 50/60Hz (CN Plug) (R290)

Cooling

Cooling of compressor 1-stage Air

Bath

Bath tank Stainless steel

Electronics

Interfaces Ethernet optional, RS232, USB

Temperature control PID1

Temperature display OLED

Temperature setting Touchscreen

Electronic Timer h:min 99 ... 59

Dimensions and volumes

Weight lbs 133.4

Barbed fittings inner diameter mm ½"

Total dimensions in. (W × L × H) 16.9 x 22 x 27.2

Body dimensions in. (W × L × H) 16.9 x 19.7 x 26.8

Filling volume l 3 ... 8

Pump connections M24x1.5 male

Temperature values

Return flow temperature max. °C +85

Working temperature range °C -20 ... +40

Temperature stability °C ±0.3

Ambient temperature °C +5 ... +40

Temperature display resolution °C 0.1

Other

Classification Classification I (NFL)

Standards CE, RoHs

Alarm notification Acoustic, Optical

IP Code IP 21

Pump function Pressure Pump

Pump type Centrifugal Pump

User Interface Language English, German

Sound pressure level dbA ² 46 ... 60

² The sound pressure level may increase as the ambient temperature rises.

Included in delivery

2 barbed fittings suitable for hoses with an inner diameter of 1/2". (Pump connections with M24x1.5 external thread)

Refer to www.julabo.com for more information regarding the entire JULABO product portfolio. Technical changes without prior notification. Images may deviate from the original. | Datasheet No.EN-US9610183/251125